


# Probe Software

Software for MicroAnalysis

Probe Software, Inc.  
885 Crest Dr  
Eugene, OR 97405 USA  
(541) 343-3400  
(541) 485-7901 (FAX)  
[sales@probesoftware.com](mailto:sales@probesoftware.com)  
[www.probesoftware.com](http://www.probesoftware.com)

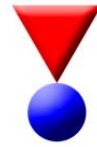


Probe for EPMA

Probe Image

## Automatic Recording of Volatile and Alternating On/Off Peak Intensities

- Automatically record all volatile calibration intensities.
- Automatically record both on and off-peak (alternating) intensity data.
- New high accuracy (quadratic exponential) fit for large volatile loss or gain.
- Display and export to ASCII file for external re-processing and plotting.

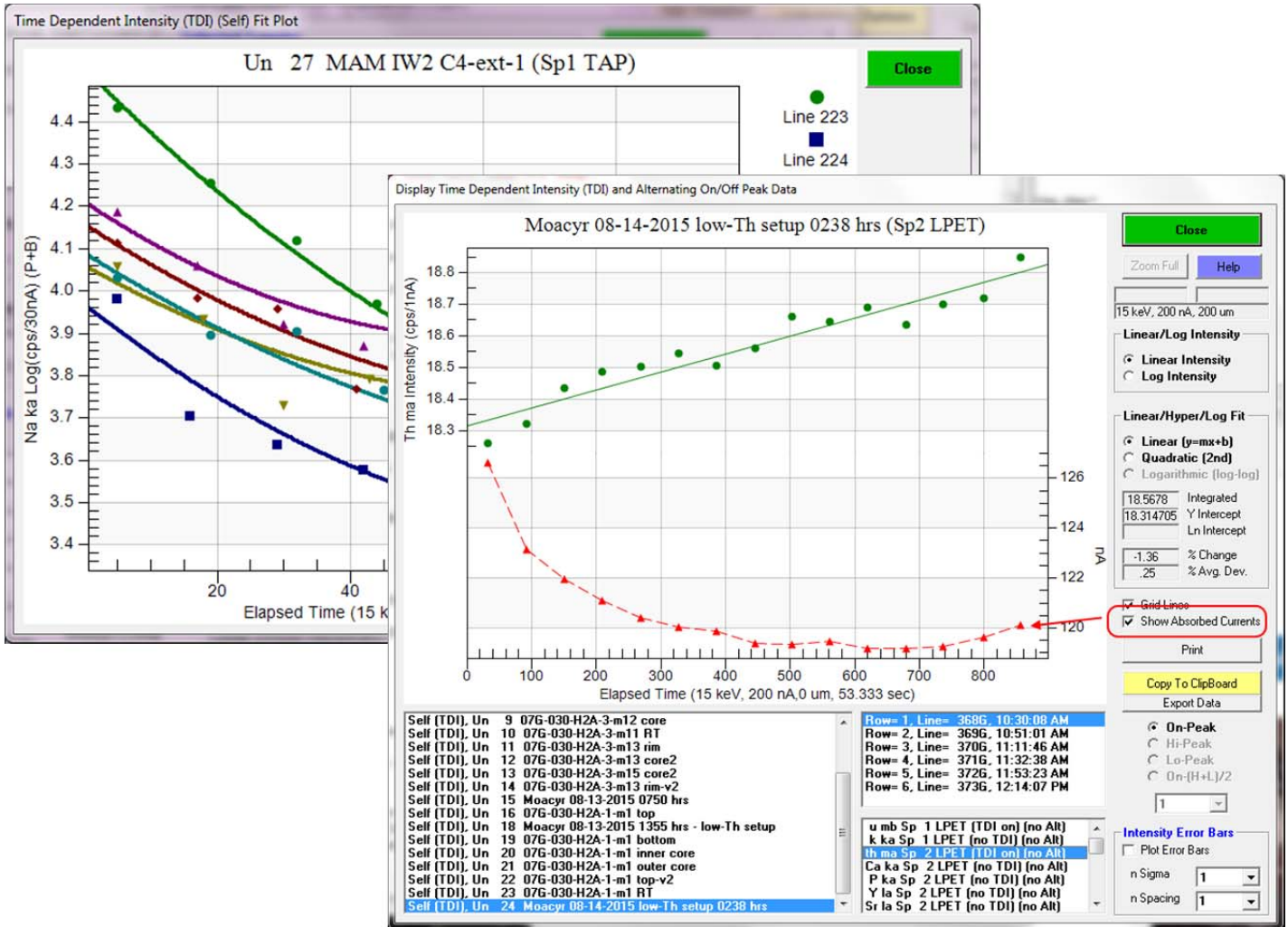

**Easily correct for sensitive samples by automatic acquisition and correction for intensity changes over time by correcting for both intensity loss or gain.**

*For extended acquisition time on trace elements, use the alternating on and off peak acquisition to integrate the peak and background intensities repeatedly and monitor long term trends in intensity changes.*