


Probe Software

Software for MicroAnalysis


Probe Software, Inc.
885 Crest Drive
Eugene, OR 97405 USA
(541) 343-3400
sales@probesoftware.com
www.probesoftware.com

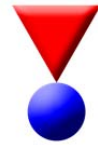

PictureSnapApp

Sample Navigation Application

- Import images from any source: slide scanner, flat bed scanner, optical microscope, SEM electron image, cell phone camera, etc.
- See your sample in true color as you navigate using the optical image as a “base map”.
- Calibrate using two reference points (or three to correct for sample rotation automatically!)
- Once your image is calibrated, double-click to drive the stage to the exact sample location.

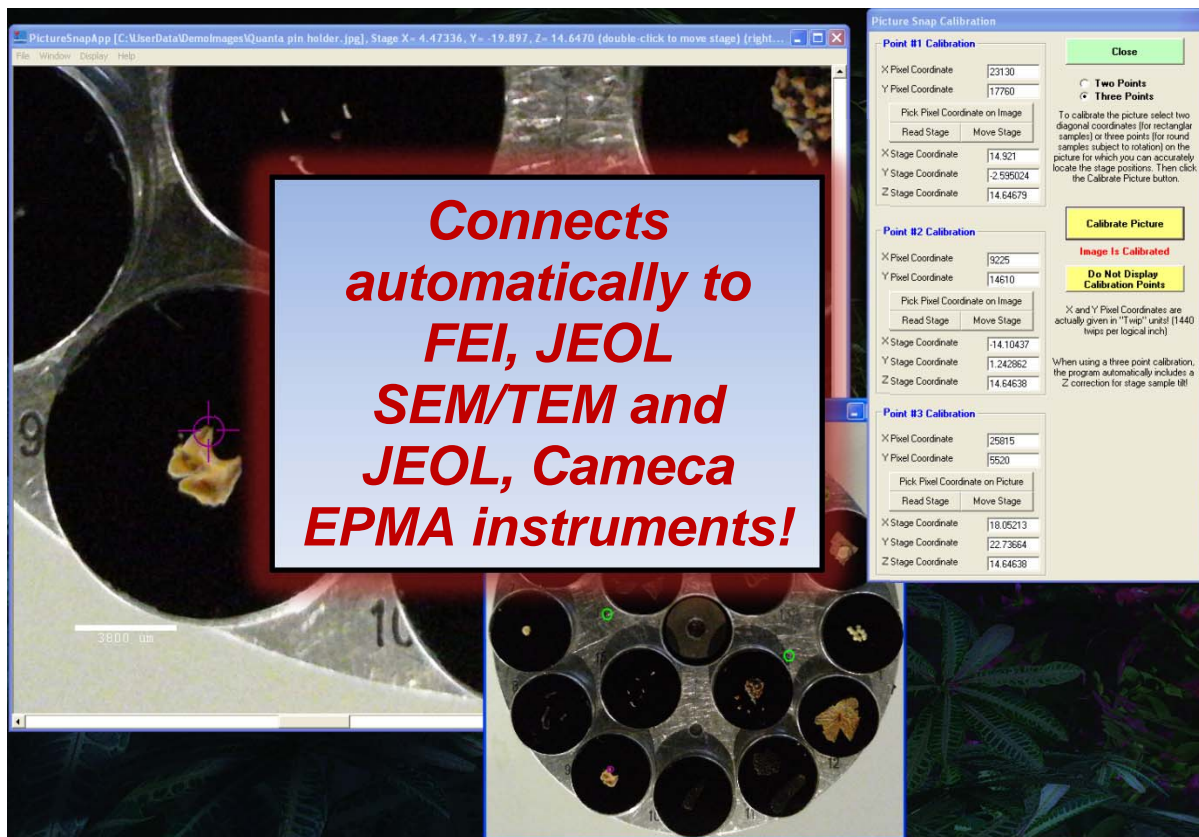

- Right-click to create unlimited sample annotations with stage coordinates that are automatically saved and exportable.
- Free download for Demonstration and Text based stage coordinate input modes!

<http://probesoftware.com/download/PictureSnapApp.msi>

- FEI/JEOL/Cameca hardware enabled PictureSnapApp can be yours for only US\$499. Please inquire about support for other instrument manufacturers.

Contact: Barbara Donovan, barbara@probesoftware.com or call (541) 343-3400 for more information.